

medisana®

HM 330

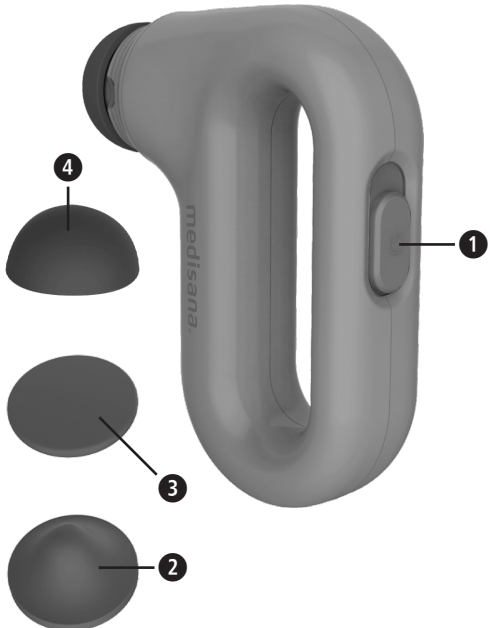Manual

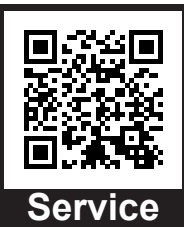Service

medisana
Carl-Schurz-Str. 2,
41460 NEUSS, Germany

REF

88261

EAN

4015588 88261 6

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Safety Instructions

Read the instruction manual carefully before using this device, especially the safety instructions, and keep the instruction manual for future use. These brief instructions are an integral part of the instructions for use. The complete instructions for use are available for download at <https://www.medisana.com>. Should you give this device to another person, it is vital that you also pass on the instructions for use.

- Please ensure that the polythene packing is kept away from the reach of children! Risk of suffocation!

- The device is intended for private use only. The appliance must not be used for commercial or medical purposes.
- If you have any health concerns, consult your doctor before use.
- Use the unit indoors only! Do not use the device in damp rooms (e.g. when bathing or showering).
- Use the device only for its intended purpose according to the instructions for use. In case of misuse, the warranty becomes void.
- If you feel pain or find the massage uncomfortable, stop the application and consult your doctor.
- You must not use the device if you suffer from one or more of the following illnesses or complaints: Circulatory disorders, varicose veins, open wounds, bruises, skin tears, phlebitis or thrombosis.
- Do not use the appliance if you are pregnant. Special care is necessary when the appliance is used on, by or near children, sick and helpless persons.
- Do not use the device near the eyes or other sensitive parts of the body.
- Do not use this device to support or replace medical applications. Chronic conditions and symptoms could worsen.
- Do not use the device while you are sleeping or in bed.
- Do not use the device before falling asleep. The massage has a stimulating effect.
- Do not use the device while driving a vehicle or operating machinery.
- This product is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the product in a safe way by a person responsible for their safety.
- This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
- Children must not play with the appliance.
- Cleaning and user maintenance must not be carried out by children without supervision.
- Do not allow the unit to come into contact with pointed or sharp objects.
- Do not cover the unit when it is switched on. Never use it under blankets or pillows. There is a risk of fire, electric shock and injury.
- Do not immerse the unit in water or other liquids.
- If liquid does enter the unit, switch it off immediately and remove the batteries. Wait until the unit is completely dry.
- Always switch off the appliance with the power switch ❶ after use.
- Do not use the appliance if any damage is visible, if it is not working properly, if it has been dropped or immersed in water. To avoid hazards, send the unit to the service centre for repair.
- Do not insert any objects into the openings of the unit.
- You may only clean the unit. In case of malfunctions, please do not repair the unit yourself, as this will automatically void all warranty claims. Contact your specialist dealer and have the repairs carried out by the authorised service.
- Only use accessories intended for the unit and supplied by the manufacturer.

SAFETY NOTES FOR BATTERIES

- Do not disassemble batteries!
- Never leave any low battery in the battery compartment since it may leak and cause damage to the unit.
- Increased risk of leakage! Avoid contact with skin, eyes and mucous membranes! If battery acid comes in contact with any of these parts, rinse the affected area with copious amounts of fresh water and seek medical attention immediately!
- If a battery has been swallowed, seek medical attention immediately!
- Replace all of the batteries simultaneously!
- Only replace with batteries of the same type, never use different types of batteries together or used batteries with new ones!
- Insert the batteries correctly, observing the polarity!
- Remove the batteries from the device if it is not going to be used for an extended period of 3 months or more!
- Keep batteries out of children's reach!
- Do not attempt to recharge batteries! There is a danger of explosion!
- Do not short circuit! There is a danger of explosion!
- Do not throw into a fire! There is a danger of explosion!

Cleaning and care

- Before cleaning the item, remove the batteries and let the unit cool down for at least 15 minutes.
- Clean the unit with a soft cloth slightly moistened with water. Never use harsh cleaning agents, petrol, thinner or alcohol.
- Then dry the unit with a dry, clean cloth.
- Do not expose the unit to direct sunlight, protect it from dirt and moisture.
- Make sure the unit is completely dry before using it again.
- Store the unit in a clean and dry place.

Voltage supply: 3 V $\approx$, 2 x 1.5V batteries (type AAA, LR03)

Disposal: This appliance must not be disposed of with household waste. Every consumer is obliged to hand in all electrical or electronic devices, regardless of whether they contain pollutants or not, to a collection point in his or her city or to a retailer so that they can be disposed of in an environmentally friendly manner. Remove the battery before disposing of the device. Do not dispose of used batteries/battery packs in household waste, but in hazardous waste or at a battery collection point in a specialist shop. Contact your local authority or dealer regarding disposal.

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