# A ILORIA

## FLASH TRAVEL

- DE** BEDIENUNGSANLEITUNG
- EN** INSTRUCTION MANUAL
- FR** MODE D'EMPLOI
- IT** ISTRUZIONI PER L'USO
- ES** MANUAL DE USUARIO
- NL** GEBRUIKSAANWIJZING
- PL** INSTRUKCJA OBSŁUGI
- LT** NAUDOJIMO INSTRUKCIJA
- CS** NÁVOD K POUŽITÍ
- DA** INSTRUKTIONSMANUAL
- FI** KÄYTTÖOPAS
- NO** BRUKSANVISNINGEN
- PT** MANUAL DE INSTRUÇÕES
- RO** MANUAL DE INSTRUCȚIUNI
- SK** NÁVOD NA POUŽITIE
- SV** BRUKSANVISNING

Version 1.0 - 10/2022

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## PRODUCT FEATURES

- Up to 40,000 brush head movements per minute
- 3 cleaning modes: Whitening - Clean - Massage
- 2 minute timer and 30 second interval for 4 quadrants
- Charging of the lithium-ion battery via USB
- Charging time approx. 8 hours
- 25 days battery life with 2x daily use of 2min
- Waterproof according to IPX6

## GETTING STARTED

Charge the toothbrush before using it for the first time. Charging time is approximately 8 hours to charge the toothbrush before first use or fully discharged battery.

Put the brush head on the toothbrush and make sure it is locked in place.

## CHARGING

If the red indicator light flashes or the product does not turn on, the battery power is insufficient and the product needs to be charged.

To charge, plug the supplied cable into the bottom of the toothbrush and connect the other side to a USB port or charging adapter.


To prevent overcharging the battery, charging will automatically stop when fully charged and the

indicator light will turn off.

Keep the charging port dry during the charging process. Close the charging port after each charging process. Make sure that the cover sits straight on the port and seals it watertight. Water ingress will permanently damage the battery.

When using for the first time or if the battery has been discharged for a longer period of time, the indicator lights up for the first 5 minutes.

## INSTALLING THE BRUSH HEAD

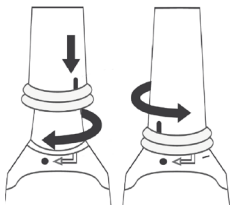Put the brush head on the neck of the toothbrush. Align the guide mark on the back of the brush head with the beginning of the guide mark on the back of the brush body.

Rotate the brush head so that the guide marks on the back of the brush head and neck are aligned.

## APPLICATION

1. Switch on the toothbrush.
2. During the entire brushing time, move the brush head slightly inclined (approx. 45°) and slowly without pressure over the teeth and chewing ridges. The brush head should be placed with one half on the gums and the other

half on the teeth to ensure cleaning even on the important transition between gums and teeth.

3. Brush your teeth for 30 seconds on each of the quadrants. The Interval Timer indicates when the quadrant should be changed by briefly interrupting. The Smart Timer ensures the two-minute brushing time recommended by the dentist.

NOTE: FLASH TRAVEL is safe to use on braces and dentures (fillings, veneers, crowns, bridges, etc.). When used on braces, the brush heads wear out more quickly.

Dentists recommend replacing brush heads every 3 months for best results.

FLASH TRAVEL has been tested and meets the safety standards for electromagnetic devices.

## **CLEANING MODES**

The FLASH TRAVEL sonic toothbrush has 3 cleaning modes: Whitening - Clean - Massage. The memory function remembers which mode you used last.

Press the power button and switch on the toothbrush. Press the power button again within 3 seconds to switch to the next mode. Repeated pressing activates all 3 modes consecutively. After the 3. and last mode, the toothbrush switches off.

The indicator lamps on the toothbrush show you which mode is activated.

<b>Mode</b>	<b>Application</b>
Clean	For daily cleaning of normal teeth and gums
Whitening	For particularly intensive cleaning of the tooth surfaces with whitening effect
Massage	For gum massage with gentle vibrations of the brush head

## **CLEANING**

- Do not clean the FLASH TRAVEL sonic toothbrush and brush heads in the dishwasher.
- Rinse the brush head and bristles with water after each use.
- Remove the brush head and rinse the neck of the toothbrush with warm water.
- Use a clean, damp cloth to wipe the entire surface of the toothbrush.
- If you do not use FLASH TRAVEL for a long period of time, remove the charging cable, clean the toothbrush and store it in a cool, dry place and avoid direct sunlight. It is recommended to recharge the toothbrush after 6 months of non-use.

## **WARNINGS**

- Consult your doctor before using FLASH TRAVEL if you have any medical problems.
- If you have had jaw or gum surgery in the last

two months, consult your doctor before using FLASH TRAVEL.

- If excessive bleeding occurs after using FLASH TRAVEL or if bleeding continues after one week of use, consult your dentist.
- Do not clean FLASH TRAVEL or the brush heads in the dishwasher.
- Do not use attachments other than those recommended by the manufacturer.
- Discontinue use if FLASH TRAVEL appears damaged in any way (brush head, handle).
- Keep the cable away from all sources of heat.
- Do not use a brush head with damaged bristles. Deformed bristles can break off when brushing. Replace the brush head every 3 months or sooner if signs of wear appear. Be aware that improper cleaning can damage the bristles. (See Cleaning)
- Use FLASH TRAVEL only for its intended use as described in this manual. This product is designed to clean teeth and gums. If you experience pain or discomfort, please stop using FLASH TRAVEL and consult your dentist or doctor as this may indicate other health problems.
- Persons (including children) with lack of experience and knowledge or reduced mental, sensory or physical capabilities should not

use this product unless they are supervised during use or unless they have been given appropriate instructions concerning the use of FLASH TRAVEL by a person responsible for their safety.

- If your toothpaste contains peroxide, baking soda or other bicarbonate, clean your brush head and treat it thoroughly with water after each use. This will prevent plastic breakage due to prolonged exposure to these chemicals.

## **TROUBLESHOOTING**

### *The toothbrush will not switch on*

You may need to charge FLASH TRAVEL.

### *Tingling sensation when brushing*

Due to the sonic vibrations, you may experience a tingling sensation when using FLASH TRAVEL for the first times. This feeling will decrease over time as you get used to the vibrations.

### *The sonic vibrations seem weak*

The bristles may be pressed too hard against your teeth. Gently place the brush head on your teeth.

## **WARRANTY & SERVICE**

Defects due to material and/or workmanship will be repaired or replaced by AILORIA under warranty if proof of purchase is provided during the warranty period. Contact the Customer Service Centre at

service@vehnsgroup.com.

Items not covered by the warranty:

- Brush heads
- Damage caused by misuse, abuse, neglect, alteration or unauthorised repair
- Normal wear and tear including scratches, discolouration or fading

## **DISPOSAL INFORMATION**

Waste electrical and electronic equipment must not be disposed of with household waste. You must therefore not simply throw away your old appliances. They must be recycled. In order to improve waste separation, the Electrical and Electronic Equipment Act therefore stipulates that all electrical and electronic equipment must be permanently marked with the corresponding symbol, the crossed-out dustbin. Free collection points for old electrical appliances are available for you to return them. You can obtain the addresses from your city or municipal administration. Batteries and accumulators that are not enclosed in the old appliance, as well as lamps that can be removed from the old appliance without being destroyed, must be removed from the old appliances before disposal. As a consumer, you are generally obliged to hand in used batteries and accumulators as well as lamps separately at the special collection points. Make sure that the batteries and accumulators are completely empty. Then carefully remove the batteries and accumulators

from the old appliance.

The following batteries are included with this item:  
1x rechargeable battery, type Lithium battery, 1040,  
3.7V,260mAh.

### Changes and errors excepted

The manufacturer accepts no liability for damage caused by incorrect, improper or incompatible use of the product.

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© VEHNS GROUP GmbH
VEHNS GROUP GmbH, Theatinerstrasse 40-42